



GIG
CYMRU
NHS
WALES

Bwrdd Iechyd Prifysgol
Aneurin Bevan
University Health Board

Dreamstation Non Invasive Ventilation (NIV) UNIT

Patient set up & user Guide



Content

Page 2/3 Dream station set up

Page 4 Starting & Stopping Treatment

Page 5 Adding Oxygen therapy

Page 6 QR Code

Dream station set up


1. Find a safe area at the side of the bed to place the unit, e.g.:- bed side table that is clear from clutter.
2. Connect the DC plug of the power supply unit to the rear of the machine and leave plugged in. (See below picture)




Plug into the main power outlet

****NB if safe to do so leave unit plugged in and in Stand-by mode at all times.***

3. Connect tubing to machine outlet (see picture below)

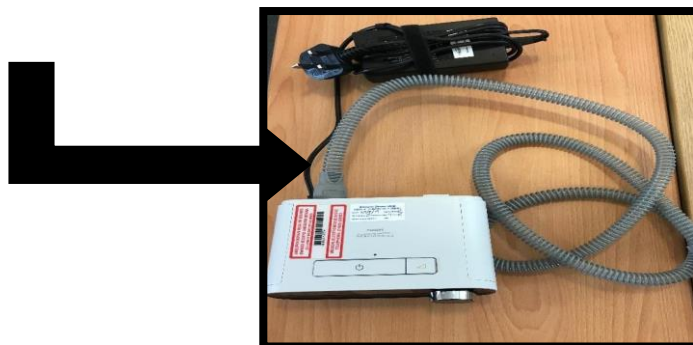

4. Connect your mask to other end of tubing (see picture)
NB Mask shown is for illustration purposes only and your own may differ from this.

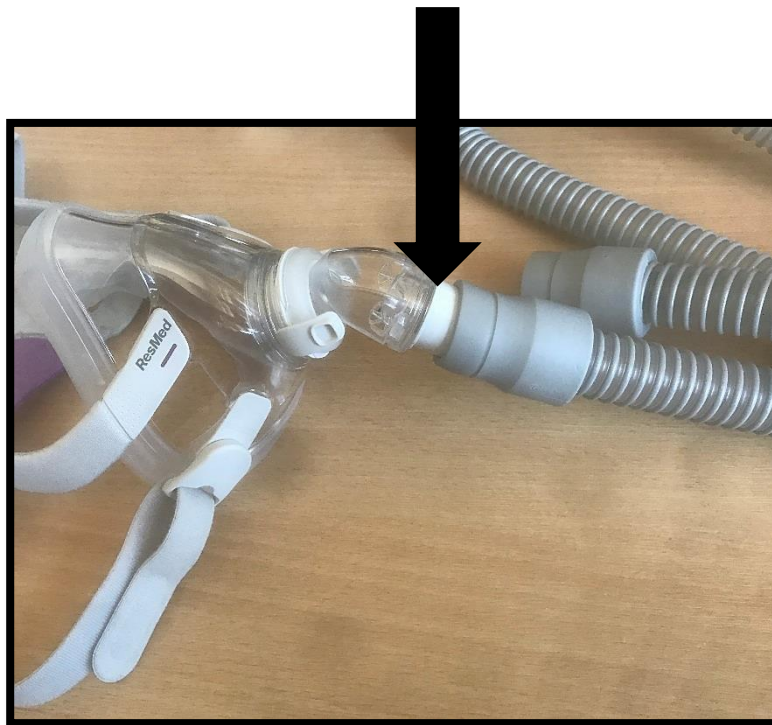

Your machine is now ready to use

Starting & Stopping Treatment

1. Ensure the mains power is on and press the start button as pictured below.
Your machine will now start to deliver air at its set pressure.



2. Apply your mask, ensuring there is a good seal, between your mask and face.
(Alternatively apply mask before starting therapy if preferred)



3. When you have completed your treatment or to stop the treatment at any time press and hold the start button as pictured above.
It may take a few seconds for the machine to stop blowing air.

Adding Oxygen Therapy



When oxygen is required you will see an additional port attached between the machine and tubing



To add oxygen remove cap as pictured



Attach oxygen tubing to port as shown

Leaflet produced by: Lyndsey Ward Respiratory ANP
Sam Jones Respiratory Nurse Specialist NIV

**"This document is available in Welsh /
Mae'r ddogfen hon ar gael yn Gymraeg".**